


POSIZIONAMENTO ELETTRODI CLICKHOLTER CARDIOLINE VETERINARIO 5 CAVI

	DESTRA	SINISTRA
SPINA SCAPOLA	Negativo BIANCO 	Negativo GIALLO 
INTERMEDIO		Terra VERDE 
STERNO	Positivo NERO 	Positivo ROSSO 

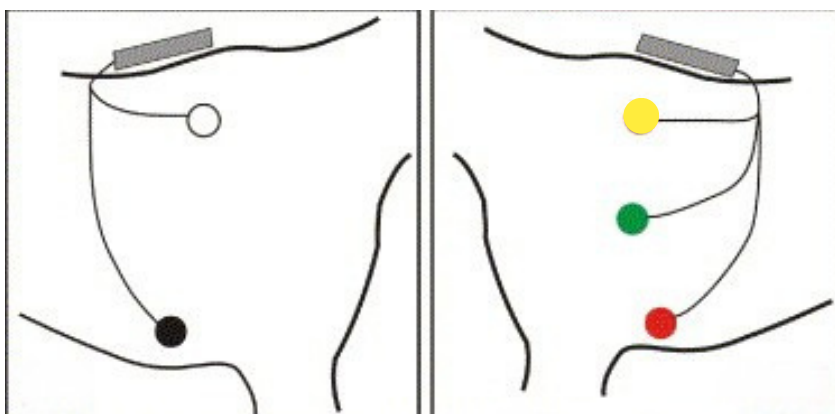

Fonte: Jean-Paul Petrie, Practical Application of Holter Monitoring in Dogs and Cats, Clinical Techniques in Small Animal Practice, Volume 20, Issue 3, 2005, Pages 173-181, ISSN 1096-2867, <https://doi.org/10.1053/j.ctsap.2005.05.006>.